# EN User Manual

## Eye Massager with Vibration Function

Model: MMOW

Thank you for purchasing our product. Before using the device for the first time, please read this manual carefully and follow all the safety instructions provided. Keep this manual for future reference, as it contains important information.

## Package Contents

1. Massage device
2. User manual

Check the product for any visible damage. If the product is damaged, it can be returned to the place of purchase for repair or replacement if the damage or defect results from a manufacturing fault. To protect the environment, dispose of the packaging in appropriate recycling containers. Keep the packaging out of reach of children. For safety reasons, do not leave any packaging components (plastic bags, cardboard boxes, polystyrene, etc.) within children's reach. Risk of suffocation!

## Intended Use

The device is designed for massaging the skin around the eyes, helping to stimulate skin nourishment and regeneration. It is intended for home use indoors only and must not be used for professional purposes. This product is intended for adult use only.

## Safety Instructions

- Read this manual carefully before using the device. Failure to follow these instructions may be dangerous or unlawful.
- Use the device only for its intended purpose.
- This product is intended for skin care around the eyes. It is not designed for laboratory, industrial, or medical use.
- The packaging is not a toy. Keep it away from children immediately after unpacking.
- The device is not a toy. Do not allow children to play with it. For adult use only.
- Repairs may only be performed by qualified manufacturer service personnel or an authorized service center. Repairs by unauthorized persons may cause damage and void the warranty.
- The manufacturer is not responsible for any damage caused by unauthorized modifications.
- Do not use the product if it is damaged.
- Protect the device from drops and strong impacts.
- Protect it from water ingress. Do not immerse in water or other liquids.
- Do not use near sinks, bathtubs, or containers filled with water.
- If an allergic skin reaction occurs during or after use, stop using the device and consult a doctor.
- Use the device at room temperature.
- Remove glasses, contact lenses, and jewelry before use.
- Wash and dry your hands before use.
- It is recommended to use the device after bathing and apply a skin care product (e.g., lotion or cream) before use.
- Do not use the device simultaneously with other similar devices.
- Avoid bringing the device too close to the eyes. Do not use on areas other than facial skin.
- The results may vary depending on skin type, usage range, and frequency.
- Do not use the device on irritated, allergic, or inflamed skin. Consult a specialist before use to determine any contraindications.
- If you notice unusual noises or smoke, immediately disconnect the power and remove the battery.
- When not in use, always protect the device with the head cover.
- After use, wipe the head with a tissue or dry cloth.
- Keep the device clean at all times. Use only a dry, soft cloth for cleaning. Do not use acidic or alkaline cleaning agents.

**WARNING:** Keep eyes closed during use. Do not apply directly to the lash line or eyelids.

**WARNING:** Do not use more than 3 times per week and no longer than 1–2 minutes per area at a time.

**WARNING:** Do not press the device too firmly against the skin or keep the head (1) in one place for too long. Move the head slowly.

## Contraindications

- Do not use on irritated skin, wounds, active acne lesions, or infections. Do not use on heavily tanned skin.
- Fever.
- Congenital heart defects or circulatory failure.
- Severe general condition or circulatory disorders.
- Pregnancy.

EN

PL

FR

DE

ES

SK

CZ

IT

RO

HU

PT

SL

GR

RU


UA

- Cancer or post-cancer surgery.
- Presence of metallic foreign bodies (e.g., dental braces, pacemakers).
- Facial nerve disorders.
- Eye irritations or infections.
- Viral, bacterial, or fungal skin diseases.
- Directly after tanning, sun exposure, or sauna use.
- After alcohol, drug use, or in a state of fatigue or physical weakness.

## Battery Safety Instructions

- The device is powered by a replaceable AA 1.5V  $\text{---}$  battery. Use only this type of battery.
- Avoid exposing batteries to extremely low or high temperatures (below 10 °C / 50 °F or above 45 °C / 113 °F). Extreme temperatures may affect battery life.
- Avoid contact between batteries and liquids or metal objects, as this may cause damage.
- Use batteries only as intended. Do not destroy, damage, or throw them into fire—this may be dangerous and cause a fire. Risk of explosion!
- Swallowing a battery can be fatal! Keep away from children and pets. Seek medical help immediately if swallowed.
- Do not disassemble, modify, open, or short-circuit batteries.
- Remove depleted batteries from the device. Leaking batteries may cause skin irritation or burns. In case of contact, wash hands with running water and consult a doctor if necessary.
- Non-rechargeable batteries must not be recharged.
- Used or damaged batteries must not be disposed of with household waste. Dispose of them at an appropriate collection point.

**NOTE:** If you are not using the device for an extended period, remove the battery.
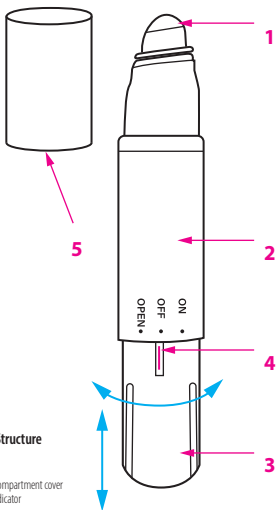
## Product Structure

1. Head
2. Handle
3. Battery compartment cover
4. Power indicator
5. Head cover

## Specifications

Power supply:	1.5 V= AA battery
Operating time:	approx. 180 min
Dimensions:	Ø 2.5 × 11.9 cm
Weight:	41 g (without battery)

## Usage

To turn the device on or off, hold the handle (2) and rotate the battery compartment cover (3) clockwise to move the power indicator (4) from OFF to ON. To turn it off, rotate the cover (3) counterclockwise from ON to OFF.

**NOTE:** It is recommended to test the device on the wrist before applying it to the face.

## Installing or Replacing the Battery

1. Hold the handle (2) and turn the battery compartment cover (3) counterclockwise until the indicator (4) aligns with the OPEN position marked on the handle (2).
2. Remove the cover, take out the old battery (if present), and insert a new one, observing polarity: positive terminal towards the handle (2), negative terminal towards the cover (3).
3. Insert the cover (3) back into the handle, aligning the indicator (4) with the OPEN mark, then turn it clockwise to the OFF position.

## Operating Instructions

1. Wash your face and ensure it is well moisturized. It is recommended to use the device after bathing.
2. Apply a suitable skincare product (e.g., lotion or cream) to the face. A facial mask is recommended for the under-eye area.
3. Turn the device on by rotating the indicator to the ON position.
4. Make sure the head (1) is not placed too close to the eyes. Do not direct it towards the eyes.
5. Note: Test the device on your wrist before use on the face.
6. Keep your eyes closed during use. Do not apply to the lash line or eyelids.
7. Hold the handle (2) with one hand and place the head (1) against the skin.
8. When used with cream or gel, distribute evenly by gently massaging from the inner to the outer eye area.
9. Keep the head moving constantly; do not keep it in one spot.
10. Massaging the cheeks upwards and outwards may help stimulate blood circulation, break down facial fat, smooth fine lines, and firm and contour the face.

## Maintenance and Storage

To clean the device, gently wipe it with a dry, soft cloth. Do not use corrosive chemicals. Do not immerse in water or any solution. Clean the head with a tissue or soft cloth. Store the device in a dry, well-ventilated place, away from heat sources. Do not store in the bathroom or humid areas. ***If not using the product for an extended period, store it in its original packaging without the battery installed.***

## Disposal

This product is marked with the crossed-out wheeled bin symbol, in accordance with European Directive 2012/19/EU on Waste Electrical and Electronic Equipment (WEEE). Products marked with this symbol must not be disposed of with household waste after the end of their service life. Users are required to take such equipment to designated collection points where it will be recycled. Proper collection and recycling help protect natural resources and have a positive impact on human health and the environment. For information on local collection points, contact your local authorities, waste management facility, or the point of sale.

According to EU Regulation 2023/1542 on batteries, this product is also marked with the crossed-out bin symbol. This means that batteries used in this product must be disposed of in accordance with local regulations and not with unsorted municipal waste. Users must use available collection systems for battery return, recycling, and disposal.

## Symbol Explanation

CE mark – The product complies with European technical standards (Directive 2014/30/EU).

Crossed-out bin – The product must not be disposed of with household waste. Non-recycled electrical equipment pollutes the air, soil, and groundwater.

Cardboard symbol – Type of packaging material: flat cardboard.

Recycling symbol – The packaging is recyclable.

Open book symbol – Read the instruction manual before use.

## Manufacturer

### TelForceOne S.A.

ul. Krakowska 119, 50-428 Wrocław, Poland

Helpline: (+48 71) 71 77 400

e-mail: info@telforceone.com; www.telforceone.pl

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